

DAICEL CHIRAL TECHNOLOGIES GLOBAL APPLICATION CENTRE

Imeglimin HCl

Column: CHIRALPAK® IK [4.6 x 250 mm, 5.0 µm]

Part No: 91325

Mobile Phase : n-Hexane: IPA: TFA: DEA (80/20/0.2/0.1, v/v/v/v)


Column Oven Temperature: 40°C

UV Detection: 240 nm

Flow Rate: 1.0 ml/min

Diluent: n-Hexane : Ethanol (1:1)

Sample Concentration: 1.0 mg/ml

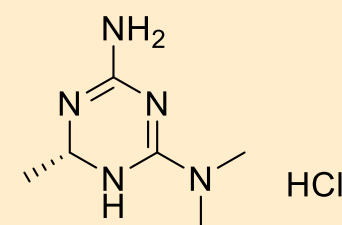

tR₁: 7.4

N₁: 4099

T₁:1.1

tR₂: 9.4

N₂:4087

T₂:1.0

Rs:3.8

